




Spectrocolorimeter specbos 4001

specbos 4001 is a very compact spectro colorimeter. It includes a small light source, a 45° 0° fiber-optic measuring head and a miniaturized spectrograph. **specbos 4001** not only allows to distinguish between different basic colors, but it can measure tiny color differences of samples. The instrument has a USB (virtual COM port) interface that can be easily installed. It needs no extra power supply.

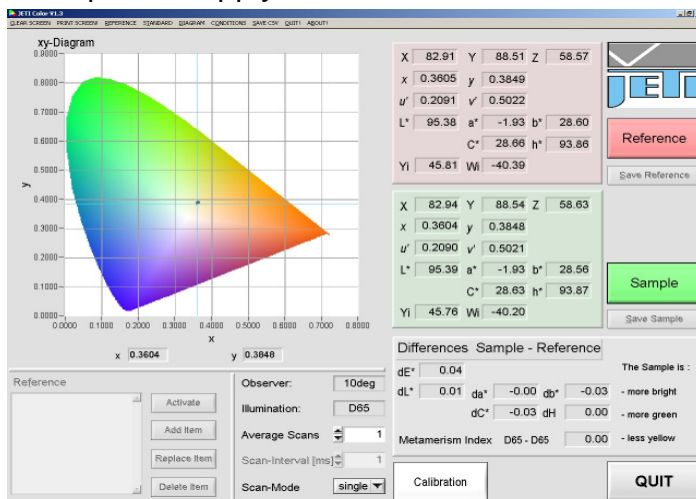

Screenshot of the JETI color software

Advantages:

- Compact instrument
- Easy installation and operation
- USB powered
- Start of measurement with external trigger

An optional holder for precise adjustment of the instrument can be delivered. It allows to align a target mark before the instrument will be pushed on the sample.

specbos 4001 can be operated with the intuitive measuring software JETI Color (for a demo version see www.jeti.com). Single and continuous measurement modes and averaging are possible. The obtained data can be exported to CSV and to ASCII files.

Specifications

Application	Surface color measurement
Optical parameters	
Spectral range	380 nm ... 780 nm or 400 nm ... 700 nm
Optical bandwidth	5 nm FWHM
Wavelengths resolution	1 nm
Digital electronic resolution	14 bit ADC
Measuring ranges and accuracies	
Measuring geometry	45°/0° (illumination 45°, measurement 0°)
Measuring probe	Diameter 18 mm
Measuring diameter	4 mm
Measuring time	typ. 0.8 s
Photometric range	0 ... 130 %
Repeatability	typ. 0.2 dE* to calib. Standard
Other technical data	
Light source	Krypton lamp with filter
Spectrometer	Imaging grating (flat field)
Light receiving element	Photodiode array 1024 pixel
Power supply	USB powered
Interface	USB 2.0 fullspeed
Dimensions	175 mm x 58 mm x 34 mm
Weight	350 g
Operating conditions	Temperature 10 ... 40 °C Humidity < 85 % relative humidity at 35 °C
Accessories (included)	PC software JETI Color for Windows 7/ XP/ Vista DLL, LabView VI's USB cable and trigger connector Operating instructions Calibration standard Calibration certificate Transport box
Calibration	BAM traceable (Bundesanstalt für Materialprüfung)
Recommended recalibration interval	1 year

Software: JETI Color

- Working under Windows 7/ XP/ Vista
- Illumination D65, A, C
- Observer 10° and 2°
- Color systems XYZ, L*a*b*(CIE), Luv, L*C*h*
- Color difference dE* (dl*, da*, db*, dC*, dh*) and metameric calculation
- Management of customer specific color references
- Averaging scans
- Continuous measurement
- Data export in CSV files
- Different color diagrams with zooming functions

JETI Technische Instrumente GmbH
Tatzendpromenade 2
D-07745 Jena

Tel. +49 (0) 3641 22 56 80
Fax. +49 (0) 3641 22 56 81
e-mail: sales@jeti.com
Internet: www.jeti.com